

 Belmont[®]

X-CALIBUR[™]

DENTAL CHAIR



Technological Innovation

Belmont has set new standards in the health field with its low-pressure electro-hydraulic devices

Over 50,000 models of X-Calibur have been sold since 1987!

The Choice is Yours, an Upholstery Design for Your Operating Position

| Models | Operating Positions | | | | |
|--------------------|---------------------|-------|-------|------|------|
| | 10:00 | 11:00 | 12:00 | 1:00 | 2:00 |
| Upholstery Design | | | | | |
| Benchmark | | X | X | X | |
| Seamed with Slings | X | X | X | X | X |



- Excellent lumbar support
- The back design allows you to get in close to the oral cavity at the 11:00, 12:00, and 1:00 o'clock positions

Benchmark: Standard Upholstery

- The beauty of UltraLeather™ from UltraFabrics
- The look and feel of calfskin

UltraLeather™ Elite Style is available at an up charge



- Slings provide arm support in the reclined positions
- Access to the oral cavity is improved at the 10:00 o'clock through 2:00 o'clock position

Seamed upholstery with slings standard



- Touch-Pad controls pre-set, auto return, and back and base movements.

Compatibility

X-Calibur Delivery Systems
for Mounting

- Over-the-patient units
- Swing-mount units
- L-R Rear mounted vacuum instrumentation




- Belmont hydraulic system
A 300 lb lifting capacity



- For proper asepsis practices. The low profile foot control activates a pre-programmed operative position, auto return, and separate base and backrest adjustment.



- Articulating dual axis prosthetic headrest adjusts to match patient's anatomy, locks securely in any desired position.

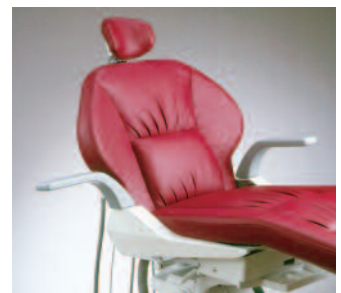

- Swingaway armrests for easy access and support on exiting the chair.

X-CALIBUR™

UPHOLSTERY AND SPECIFICATIONS

ELITE UPHOLSTERY:

UltraLeather™ provides the look and feel of real leather combined with durability and abrasion resistance.

UltraLeather™ is a trademark of UltraFabrics, LLC

BENCHMARK AND SEAMED WITH SLINGS UPHOLSTERY


Duration vinyl offers —

- Outstanding clean-ability.
- Chemical resistance.
- Bacterial resistance.

The Mor-gard™ formulation from Morbern.

It is the responsibility of the user to verify that the manufacturer of the disinfectant states that their product is safe to use on vinyl and plastic materials.

SPECIFICATIONS:

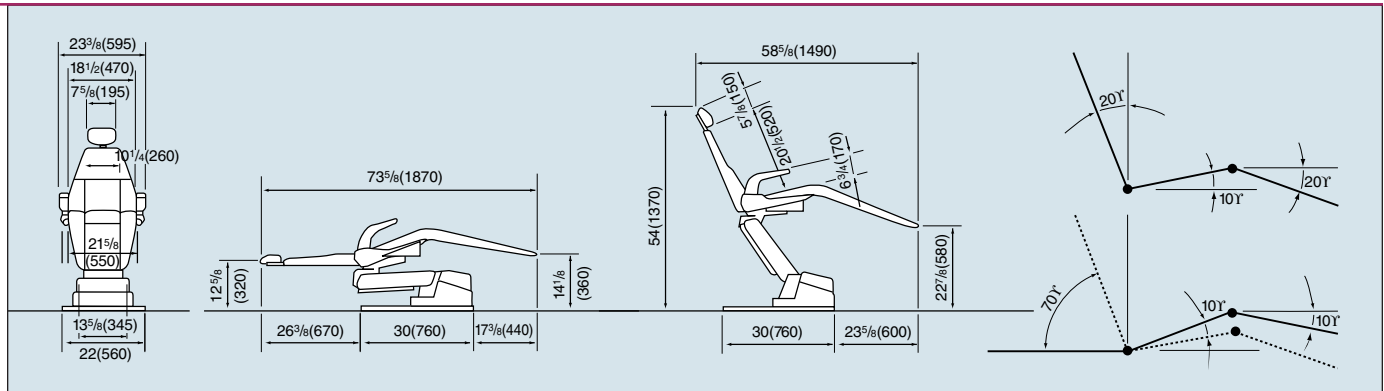

- Power sealed hydraulic system powered by 4.6A motor pump
- Base plate: 1/2" steel
- Bearings at link points — low friction insteel jacket
- Steel seat cradle
- Steel backrest support
- Electrical requirements: 114 VAC, 4.6A

- Hospital grade plug
- Base component; housing formed from solid color acrylonitrile butadiene styrene (ABS)
- All exposed ferrous parts covered with corrosion-resistant paint or plating
- Weight: 381 lbs. (173 kg)

Warranty: 5 Years on chair
1 Year on upholstery



ETL listed
Models, specifications and prices subject to change without notice
BRH00006 REV010808

Belmont Equipment
A Division of Takara Belmont USA, Inc.
101 Belmont Drive, Somerset, NJ 08873
(732) 469-5000 FAX: (800) 280-7504

1-800-223-1192
www.belmontequip.com