

Swing-Mounted Delivery With Optional Assistant's Vac Pac

- BDS2575 Delivery System
- BDS0050 Assistant's Vac Pac (For B50N, B50 & B20N Chairs.) **Doctor's Control**
- 3 Handpiece automatic control Individual handpiece pressure adjustment Individual handpiece water adjustment Master handpiece air coolant adjustment Non-retracting water coolant with master ON/OFF Handpiece flush system 4-hole Midwest style straight gray tubing 3-way syringe Wet/Dry foot control Stainless steel tray with non-slip pad

Arm System

Balance arm with air brake Delivery swing arm with city/bottled water system

Chair Adapter

Chair adapter for B50N, B50 or B20N Adapter for Quolis® chair (upcharge option)

Junction Box

J-Box with air and water floor utilities 5' Umbilical

Assistant Instruments (option)

Rear-mounted swing arm 4-position crescent shaped instrument holder 3-way syringe, HVE and saliva ejector 1/4'' water QD with outlet flow control Solids collection canister with disposable screen

Shown with BDS0050 assistant's Vac Pac option on a B50 chair

Swing-Mounted Delivery With Optional Monitor Mount & Optional Assistant Vac Pac

- BDS2575 Delivery System
- BDS2589 Monitor Arm & FPMM
- BDS0050 Assistant's Vac Pac (For B50N, B50 & B20N Chairs.)

Doctor's Control

3 Handpiece automatic control Individual handpiece pressure adjustment Individual handpiece water adjustment Master handpiece air coolant adjustment Non-retracting water coolant with master ON/OFF Handpiece flush system 4-hole Midwest style straight gray tubing 3-way syringe Wet/Dry foot control Stainless steel tray with non-slip pad

Arm System

Balance arm with air brake Delivery swing arm with city/bottled water system Monitor swing arm (BDS2589 & FPMM options)

Chair Adapter

Chair adapter for B50N, B50 or B20N Adapter for Quolis® chair (upcharge option)

Junction Box

I-Box with air and water floor utilities 5' Umbilical

Assistant Instrument Arm (option)

Rear-mounted swing arm 4-position crescent shaped instrument holder 3-way syringe, HVE and saliva ejector 1/4'' water QD with outlet flow control Solids collection canister with disposable screen



Shown with optional accessories : Unit mounted touchpad, BDS0050 assistant's Vac Pac, BDS2589 monitor arm & FPMM flat panel monitor mount on a B50 chair (Monitor not available)


CONFIGURATION AND MODEL NUMBER

	ONFIGUR	ATION AND MODEL NOMBER				
			Over-the-Patient Delivery		Cuspidor (only) w/ Asst. Instr.	
Components			BDS2500 BDS2510		BDS2520	
Doctor's	control					
Post mou	unt utility sect	ion (PMU) with city/bottle water system				
Gravity	drain cuspido	r				
Telescop	oing assistant	instrument arm with syringe, HVE & SE				
Dental light			(*Note)	(*Note)	(**Note)	
(*Note)	Choose from	n the following:				
	BDS2536N	Bel-Halo LED light with curved light post	AL901R	Bel-Halo L	ED light and B l	
	BDS2537N	Clesta LED light with curved post	AL301R	Clesta LEE	light and BDS	
	BDS2535N	Clesta halogen light with curved light post	AL511T	gen light and B		
(**Note)	FHAA03A0	Light post base assembly required for mount	ing all lights o	nto BDS2520	cuspidor	

(**Note) FHAA03A0 Light post base assembly required for mounting all lights onto BDS2520 cuspidor Swing-Mounted Delivery					Monitor Arm (For BDS2575)	Assistant's Vac Pac For B50N, B50 & B20N Chairs
Components	BDS2570	BDS2572	BDS2573	BDS2575	BDS2589	BDS0050
Doctor's control						
Left / right swing arm with city/bottle water system						
Left/right curved post light	Clesta halogen	Bel-Halo LED	Clesta LED			
Left/right rear mounted assistant instrument arm						
Left / right monitor swing arm					***	

(***Note) Order BDS2589 monitor swing arm and FPMM monitor mount for attaching to BDS2575 Delivery only. For stand-alone monitor arm-without delivery, choose **EVG3590B** (includes pivot assembly, chair adapter flat panel monitor mount.)

DIMENSIONS

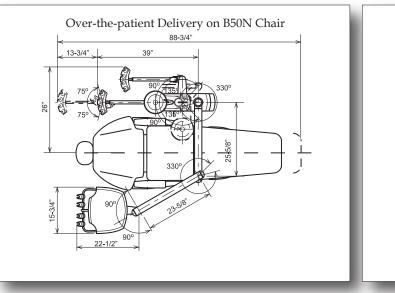

Models, specifications, prices and warranty policies subject to change without notice BRH00052 REV021419


Takara Company Canada Ltd. 2076 South Sheridan Way Mississauga, ON L5J 2M4, Canada (905) 822-2755 (800) 268-5351 dental.takarabelmont.ca

BDS2575 Delivery System




Belmont

DS0012 24" light post **0012** 24" light post **DS0012** 24" light post

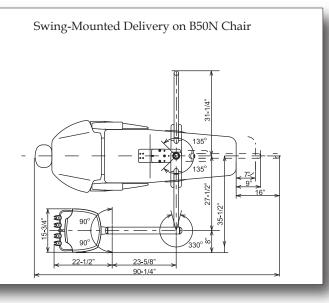

Belmont

Belmont Equipment Division of Takara Belmont USA, Inc. 101 Belmont Drive, Somerset, NJ 08873 Phone : 1-800-223-1192 www.belmontequip.com

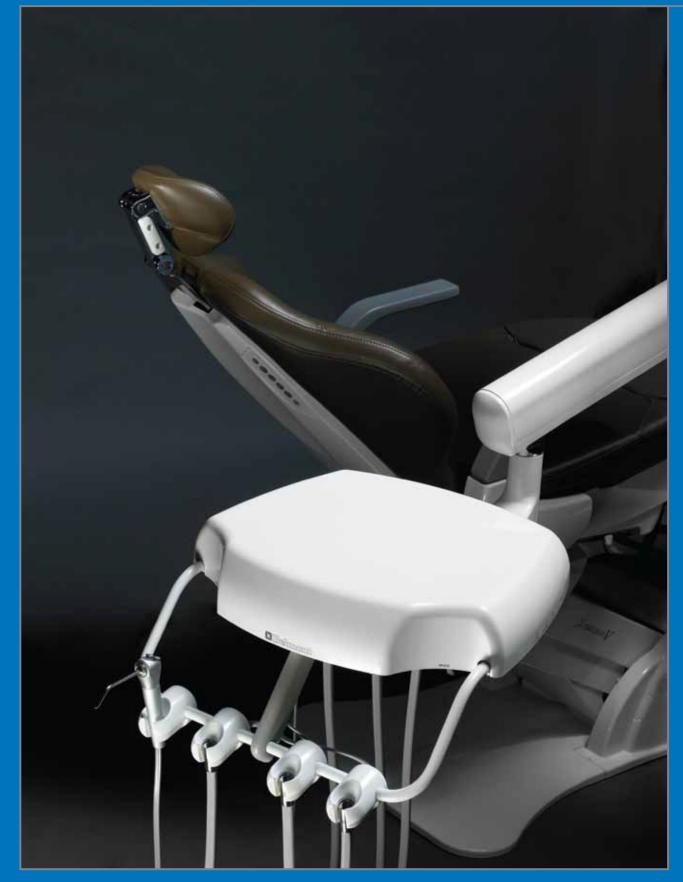

X-CALIBUR **Delivery Systems**

Over-the-Patient Delivery Configurations

X-Calibur over-the-patient delivery systems are available in 2 separate configurations to meet the needs of traditional dentistry and beyond

Over-the-Patient Delivery, With Cuspidor

• BDS2500 With Assistant Instruments

Doctor's Control

3 Handpiece automatic control Individual handpiece pressure adjustment Individual handpiece water adjustment Master handpiece air coolant adjustment Non-retracting water coolant with master ON/OFF Handpiece flush system 4-hole Midwest style straight gray tubing 3-way syringe Wet/Dry foot control Stainless steel tray with non-slip pad Remote toggles for cuspidor cupfill & bowl rinse

Arm System

Balance arm with air brake Horizontal arm

Post Mounted Utility with Cuspidor

Post mounted utility (PMU) side box City/bottled water system Swing-away cuspidor upper section Removable glass bowl for easy cleaning Push buttons for cupfill and bowl rinse Adjustable timed bowl rinse 1/4" Water QD with outlet flow control 1/4" Air OD Solids collection canister with disposable screen

Assistant Instruments

Telescoping assistant instrument arm 4-position crescent shaped instrument holder 3-way syringe, HVE and saliva ejector

Chair Adapter

Chair adapter for B50N, B50 or B20N Adapter for Quolis® chair (upcharge option)

Junction Box

5' umbilical I-Box with air and water floor utilities

Light (option)

Choose Bel-Halo LED, Clesta LED or Clesta halogen lights



FLEXIBILITY



ERGONOMICS

Over-the-Patient Delivery, No Cuspidor

• BDS2510 With Assistant Instruments

Doctor's Control

3 Handpiece automatic control Individual handpiece pressure adjustment Individual handpiece water adjustment Master handpiece air coolant adjustment Non-retracting water coolant with master ON/OFF Handpiece flush system 4-hole Midwest style straight gray tubing 3-way syringe Wet/Dry foot control Stainless steel tray with non-slip pad

Arm System

Balance arm with air brake Horizontal arm

Post Mounted Utility

Post mounted utility (PMU) side box City/bottled water system 1/4'' Water QD with outlet flow control 1/4" Air QD Solids collection canister with disposable screen

Assistant Instruments

Telescoping assistant instrument arm 4-position crescent shaped instrument holder 3-way syringe, HVE and saliva ejector

Chair Adapter

Chair adapter for B50N, B50 or B20N Adapter for Quolis® chair (upcharge option)

5' umbilical

Light (option) Choose Bel-Halo LED, Clesta LED or Clesta halogen lights

BDS2500 Delivery System with Cuspidor & Telescoping Assistant Instrument Arm Shown with touchpad option, BDS2536N curved post version of Bel-Halo LED light, mounted on a B20⁺chair **Junction Box** I-Box with air and water floor utilities

Belmont

RELIABILIT



BDS2510 Delivery System with Telescoping Assistant Instrument Arm (No Cuspidor) Shown with AL511T Clesta halogen light option and BDS0012 24" light post option, mounted on a B20⁺chair

Swing-Mounted Delivery Configurations

X-Calibur swing-mounted delivery systems meet the needs of practices that require flexibility to reposition operator set-up from the right side of the chair to the left side within seconds

Swing-Mounted Delivery, With Dental Light and Optional Assistant Vac Pac

- BDS2570 Delivery with Clesta halogen Light
- (BDS2572 Delivery with Bel-Halo LED Light, or BDS2573 Delivery with Clesta LED Light)
- BDS0050 Assistant's Vac Pac (For B50N, B50 & B20N Chairs)

Doctor's Control

3 Handpiece automatic control Individual handpiece pressure adjustment Individual handpiece water adjustment Master handpiece air coolant adjustment Non-retracting water coolant with master ON/OFF Handpiece flush system 4-hole Midwest style straight gray tubing 3-way syringe Wet/Dry foot control Stainless steel tray with non-slip pad

Arm System

Balance arm with air brake Delivery swing arm with city/bottled water system Light swing arm assembly

Light Choices

Bel-Halo LED light (2970 FTC - 370 FTC Variable) 5000°K approximates natural daylight at <u>all</u> settings Clesta LED light (2601 FTC - 288 FTC Variable) 5000°K approximates natural daylight at all settings Clesta light (2600 FTC, 2000 FTC and 740 FTC)

5000°K approximates natural daylight at high setting

Chair Adapter

Chair adapter for B50N, B50 or B20N Adapter for Quolis® chair (upcharge option)

Junction Box

J-Box with air and water floor utilities 5' Umbilical

Assistant Instruments (option)

Rear-mounted swing arm 4-position crescent shaped instrument holder 3-way syringe, HVE and saliva ejector 1/4" water QD with outlet flow control Solids collection canister with disposable screen



BDS2570 Swing-Mounted Delivery System with Swing-Mounted Clesta halogen Light Shown with BDS0050 rear-mounted assistant's Vac Pac & 90175 additional HVE options, mounted on a B50 chair