



Clesta Halogen Dental Light

■ Clesta Halogen

Dental Light

A Belmont staple in delivering simplified operations and efficiency.

Easy Positioning

Convenient ergonomically shaped handles permit easy positioning of the light head for directing illumination where it is needed.

Asepsis Potential

Touchless sensor facilitates hands-free activation for minimizing the risk of cross-contamination of the light head surfaces.



Intensity Selection

Three-way intensity switch adjusts illumination to high intensity, medium intensity, and composite-safe mode.

Color Rendering

Using the high intensity setting, the 5000°K color temperature approximates natural daylight.

Reflective Design


Special Dichroic reflector filters out infrared rays and dissipates heat, so cool white light illuminates the oral cavity.

■ Low Effort Maintenance

Quarter-turn screw releases the front lens cover assembly for fast and easy bulb replacement.

■ Productivity Through Obstruction

Shadowless illumination technology allows for uninterrupted lighting to the patient oral cavity, even with a 50% obstruction between the light head and the patient.

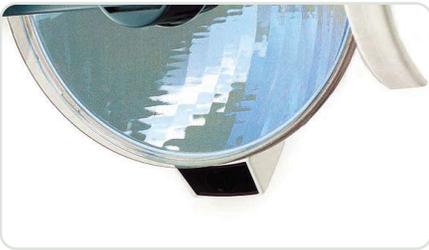

■ Mounting Options

Configurations include ceiling mount, track mount, wall mount, post mount, unit mount, swing mount and cabinet mount.

AL512T Ceiling Mount Light
(for operatories with 8 to 10 ft ceilings)



AL511T Post Mount Light
(post not included)




AL515T Track Mount Light
(for operatories with 8 to 10 ft ceilings)




BDS2530 Left - Right Swing
Mount Light




Unit Mount

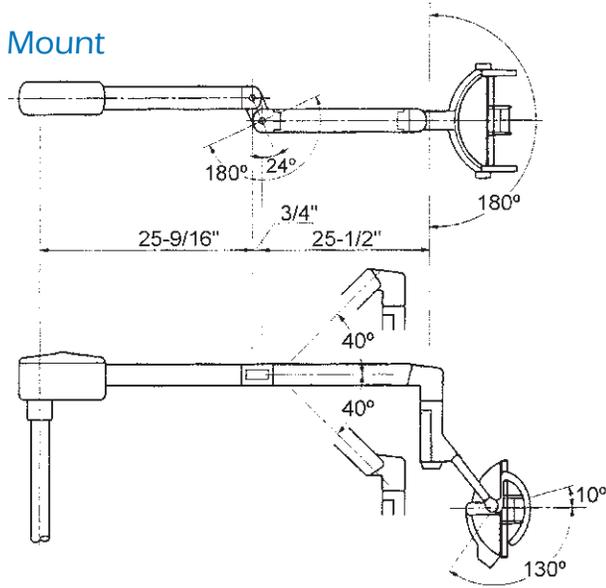


Wall Mount

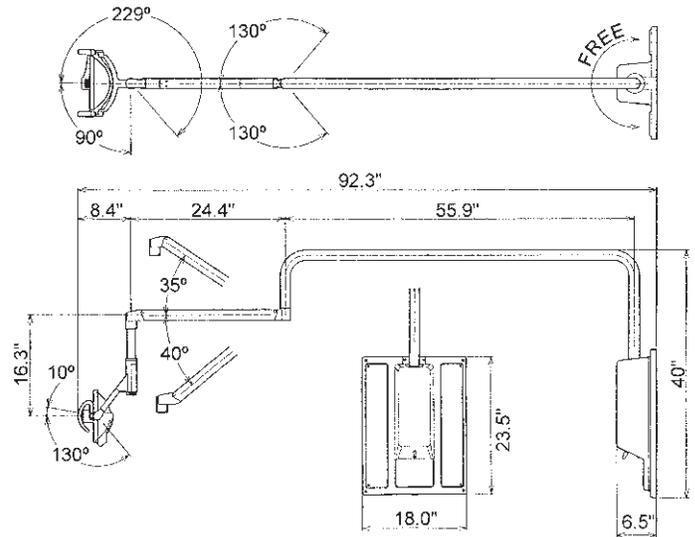

Ceiling Mount

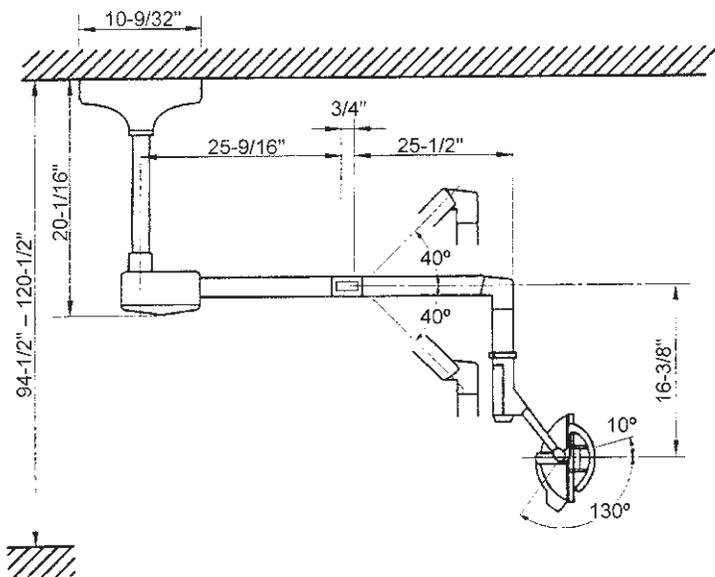

Track Mount

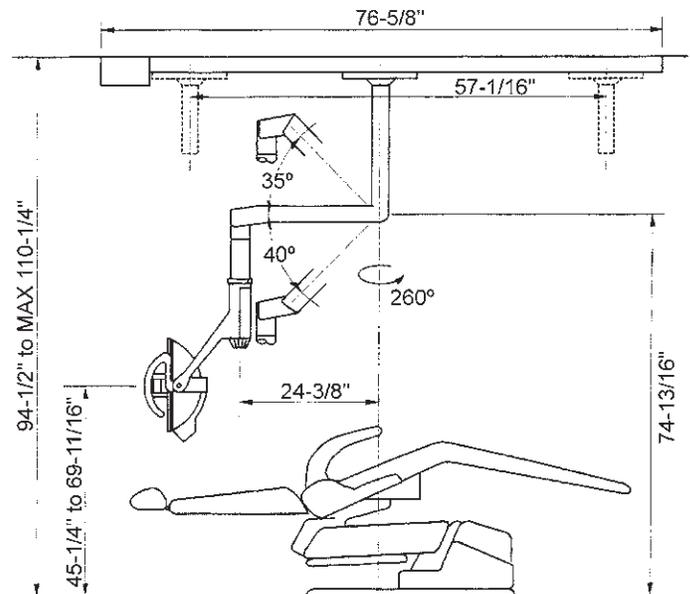

SPECIFICATIONS

Light Source:	Halogen Bulb
Focal Distance:	27-1/2" (698 mm)
Intensity:	740 - 2,600 FTC (Stepless)
Composite Safe Mode:	740 FTC
Color Temperature:	3,600° - 5,000° K
Light Pattern:	2.8" x 7.5" (71 mm x 190 mm)
Operation:	Touchless sensor
Bulb Type:	24V / 60W Tungsten Halogen
Spare Bulb:	1 Spare Bulb Included

Belmont Equipment
101 Belmont Drive
Somerset, NJ 08873
(800) 223-1192
www.belmontequip.com

Takara Company Canada Ltd.
2076 South Sheridan Way
Mississauga, ON
L5J 2M4, Canada
(800) 268-5351
dental.takarabelmont.ca

Belmont
Delivering Value in Quality

@BelmontEquipUS

