


AL108

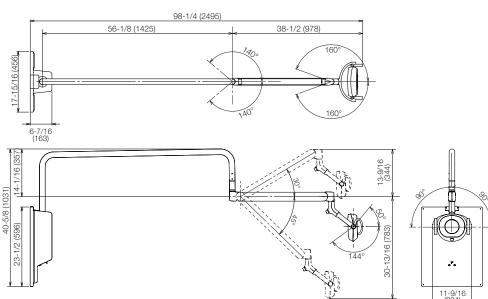

# **SPECIFICATIONS**

Light Source:	Six LEDs
In Treatment Mode	
Central Illuminance:	3,000-28,000 lx
Correlated Color Temperature:	5,000 K
In Composite Safe Mode	Compliant with ISO 9680: 2014 5.2.10
Central Illuminance:	9,000 lx
Correlated Color Temperature:	2,700 K
Standard Irradiation Distance:	25.59" (650 mm)
Pattern Dimensions:	3-3/8" x 6-1/8" (85 x 155 mm)
Operation:	Touchless sensor
Intensity Adjustment:	Stepless rotating knob
Head Positioning:	Rotation of head at three axis

**Belmont Equipment** 101 Belmont Drive Somerset, NJ 08873 (800) 223-1192 www.belmontequip.com













Shine Brighter, Work Better, Last Longer.

(3)Belmont





# Bel-Nova

# Dental Light

See your practice in a new light.

#### The Belmont Standard

Belmont's Bel-Nova is bringing you the best dental lighting technology available. This world class LED light offers better illumination in a high performance oval lighting pattern while letting off less heat and sound than traditional halogen bulbs.

#### Out of the Shadows

This cutting edge, six bulb pattern focuses more light into the patient's oral cavity, giving doctors a better vision of their patient's mouth by producing a oval pattern. It also keeps light from reflecting painfully into the patient's eyes during procedures.





### Proof In Pictures

The lighting disbursement pictured below creates less shadows, allowing doctors to see further back into their patients' mouths, leading to enhanced performance and more accurate diagnoses.





## Touchless

Keep your operatory sanitary and easy to maintain with the Bel-Nova's touchless, soft on/off sensors. Avoid germs, residue, and fatigue by simply waving your hand over the sensors to adjust the lighting, and even to switch to composite mode.

### Stepless

The days of catering to predetermined lighting settings are over. Adjust luminosity to your exact preference with the Bel-Nova's stepless technology.

#### Seamless

The Bel-Nova is made to seamlessly integrate into your existing operatory and daily procedures. With multiple set-up options featuring wider angles, longer arm reach, and 3-axis rotation the Bel-Nova and Belmont have you covered - from every position.

### Mounting Options

Unit mount (Quolis,Evogue) Color: Light Grey and White



Track mount Color: Light Grey



Post mount Color: Light Grey and White



Cabinet mount Color: Light Grey



Ceiling mount
Color: Light Grey



Wall mount Color: Light Grey



Belmont Equipment

Belmont Equipment

Belmont Equipment