


4-PILOT HOLES DRILL DIA : $\phi 2.4$ (3/32)

Note
 12AWG wire must be used
 up to a distance of 50'

CONNECT WIRES IN JUNCTION BOX AS SHOWN BELOW
 X-RAY INSTALLATION WIRING DIAGRAM

	100/110/120V	220/230/240V
L	BLACK	BROWN
N	WHITE	BLUE
\perp	GREEN	GREEN/YELLOW

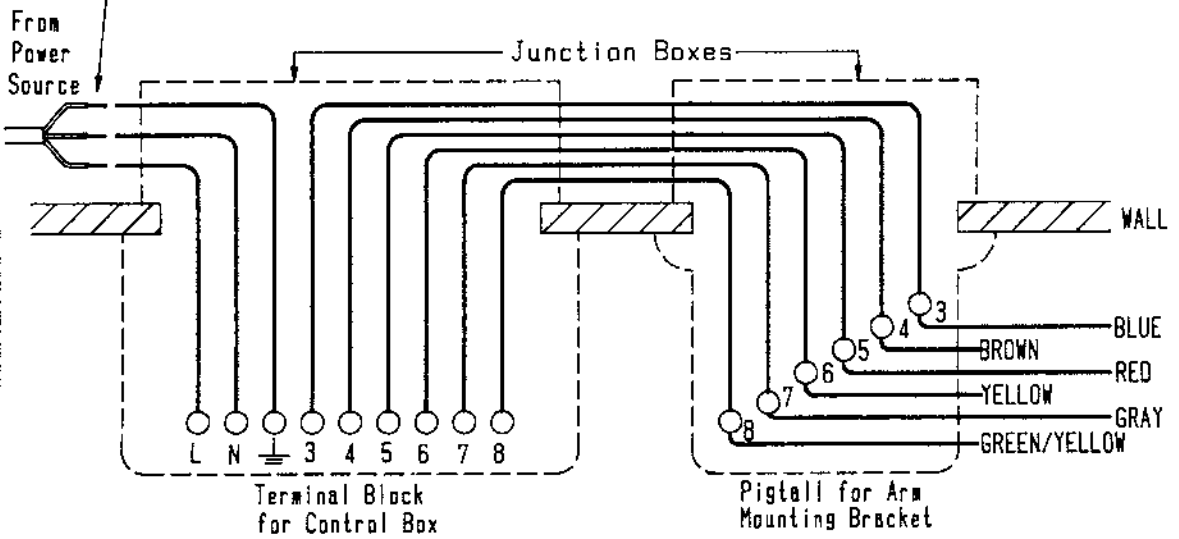


285
 (11-7/32)

13 (33/64) 130 (5-1/8) 13 (33/64)

32
 (1-1/4)

Location for Wires

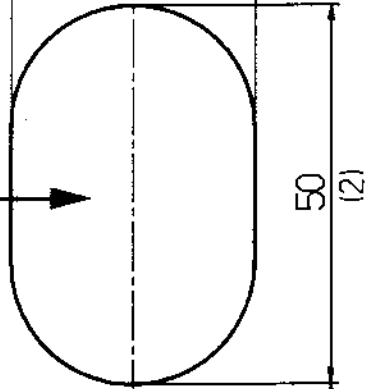

12
 (15/32)

16.2
 (5/8)

1270 (50) FROM FLOOR