




Shine Brighter,
Work Better,
Last Longer.


Bel-Nova
LED Dental Light



AL-107

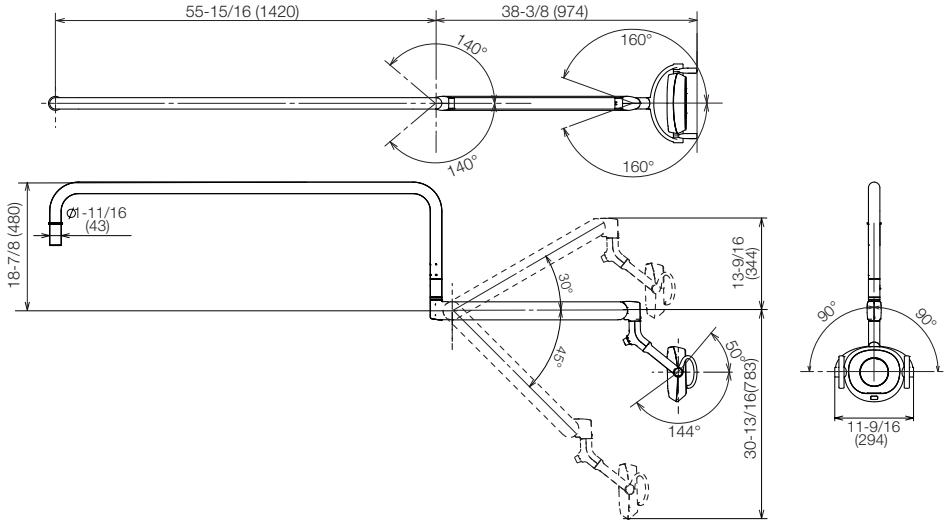

AL108

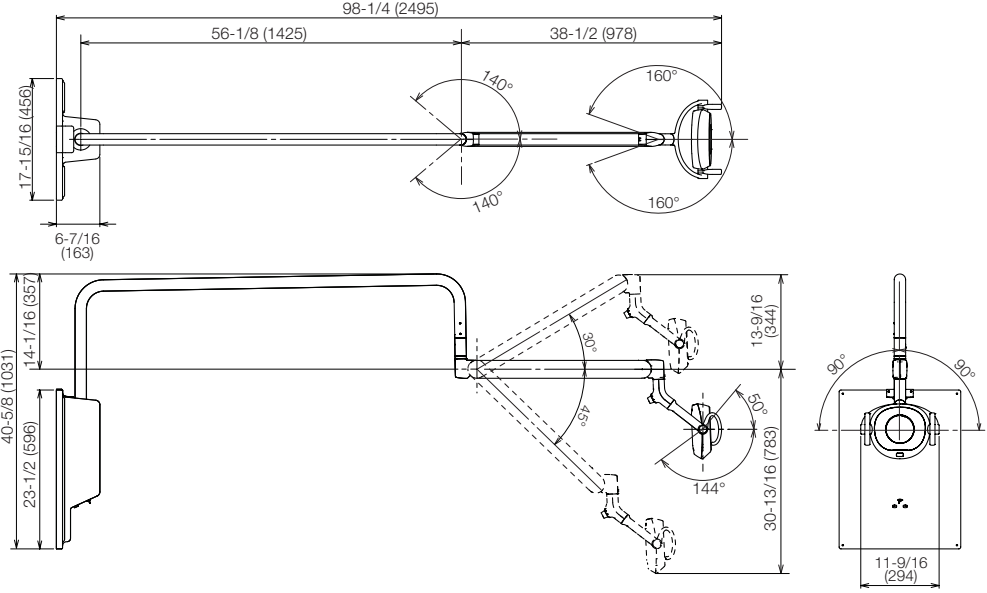

SPECIFICATIONS

Light Source:	Six LEDs
In Treatment Mode	
Central Illuminance:	3,000-28,000 lx
Correlated Color Temperature:	5,000 K
In Composite Safe Mode	Compliant with ISO 9680: 2014 5.2.10
Central Illuminance:	9,000 lx
Correlated Color Temperature:	2,700 K
Standard Irradiation Distance:	25.59" (650 mm)
Pattern Dimensions:	3-3/8" x 6-1/8" (85 x 155 mm)
Operation:	Touchless sensor
Intensity Adjustment:	Stepless rotating knob
Head Positioning:	Rotation of head at three axis

Belmont Equipment
101 Belmont Drive
Somerset, NJ 08873
(800) 223-1192
www.belmontequip.com



■ Bel-Nova

Dental Light


See your practice in a new light.

The Belmont Standard

Belmont's Bel-Nova is bringing you the best dental lighting technology available. This world class LED light offers better illumination in a high performance oval lighting pattern while letting off less heat and sound than traditional halogen bulbs.

Out of the Shadows

This cutting edge, six bulb pattern focuses more light into the patient's oral cavity, giving doctors a better vision of their patient's mouth by producing a oval pattern. It also keeps light from reflecting painfully into the patient's eyes during procedures.

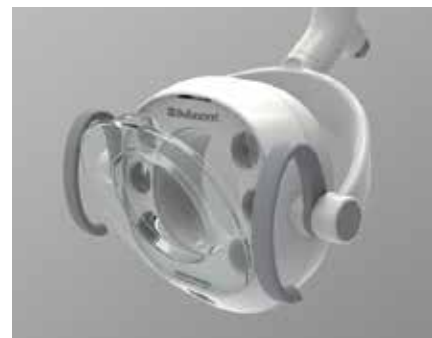

■ Proof In Pictures

The lighting disbursement pictured below creates less shadows, allowing doctors to see further back into their patients' mouths, leading to enhanced performance and more accurate diagnoses.



Less is More



■ Touchless

Keep your operatory sanitary and easy to maintain with the Bel-Nova's touchless, soft on/off sensors. Avoid germs, residue, and fatigue by simply waving your hand over the sensors to adjust the lighting, and even to switch to composite mode.

■ Stepless

The days of catering to predetermined lighting settings are over. Adjust luminosity to your exact preference with the Bel-Nova's stepless technology.

■ Seamless

The Bel-Nova is made to seamlessly integrate into your existing operatory and daily procedures. With multiple set-up options featuring wider angles, longer arm reach, and 3-axis rotation the Bel-Nova and Belmont have you covered - from every position.

■ Mounting Options

Unit mount (Quolis, Evogue)
Color: Light Grey and White



Post mount
Color: Light Grey and White



Ceiling mount
Color: Light Grey



Track mount
Color: Light Grey



Cabinet mount
Color: Light Grey



Wall mount
Color: Light Grey

